


LUMA MASK PROFESSIONAL

User Manual

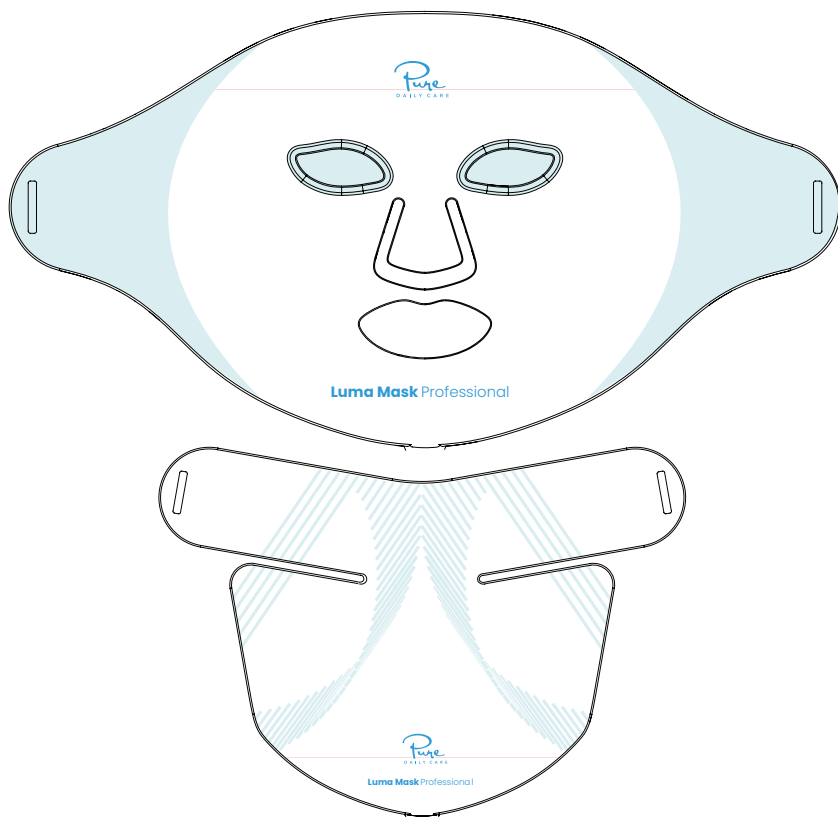

SAFETY

For your safety, please read through all warnings and additional contents listed in this manual. As for all skin care devices and therapies, always consult your doctor prior to use.

Avoid staring directly at the side of the mask that produces the light (the side that goes against your skin).

Keep your eyes closed during time of use to protect your eyes from unwanted light exposure.

Use the mask(s) inside only to reduce exposure to the sunlight.

Keep the mask(s) away from water. Do not use while bathing or where you could accidentally drop the mask into water.

Always consult a physician if you are pregnant, nursing, under the age of 21 years or have skin concerns of a medical nature.

It is recommended that you remove all jewelry before using the Luma Mask Professional.

If you experience pain or discomfort during your session, immediately stop using the Luma Mask Professional and consult your physician.

MAINTENANCE

To keep your Luma Mask Professional clean, always wipe with a damp cloth after each use. For an occasional deep clean, gently wipe masks with a small amount of 70% alcohol. Store flat to maintain the longevity of the masks.

Part Information

Face part size: 39x21x5.3cm
Neck part size: 22x30x5.3cm
Product weight: 535g
Package size: 42x29.5x7cm
Package weight: 1.45g
Box size: 44.2x44x24cm

Specifications

Material: Food grade silicone
Battery capacity: 1000mAh polymer lithium battery
Input: AC100–240V 50/60HZ
Output: DC5V 1A
Power: 5w
Lamp beads for face: 105pcs
Lamp beads for neck: 114pcs
5 levels of intensity
Charging interface: Micro USB
Charging time: 40 minutes
Available time: 120 minutes
Color: White



GET TO KNOW YOUR LUMA MASK PRO

WHAT'S INCLUDED:

1. Silicone Face Mask with seven light therapy functions:


- a) Red Light: Anti-Aging
- b) Blue Light: Anti-Blemish
- c) Green Light: Reduces Pigmentation
- d) Yellow Light: Reduces Redness
- e) Purple Light: Cell Rejuvenation
- f) Cyan Light: Reduces Irritation
- g) White Light: Boosts Skincare Treatments

2. Silicone Neck Mask (Note: all seven light functions can be used across both face and neck)

3. Adjustable straps

4. Eye area comfort inserts (2 sets)

5. Mask light controller

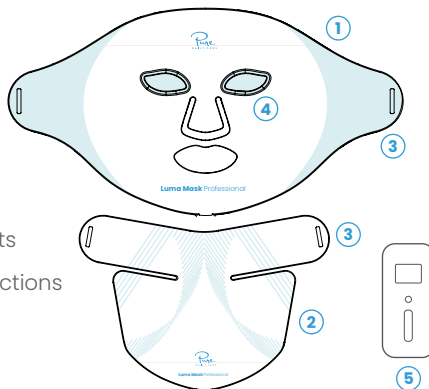

HOW TO USE YOUR LUMA MASK PRO

1. Remove the Luma Mask Professional parts from packaging and lay on a flat surface.

2. Start by connecting the supplied USB cord to the controller.

3. To secure the masks to your face and neck, use the adjustable straps to comfortably fit both masks to your face and neck.

4. To begin the light therapy, select your light of choice, followed by the length of time.

Note: It's recommended that you use it for 15-20 minutes up to 3 times per week.

5. You can adjust the strength of the light by pressing power button to adjust to your preference.

6. Use the mask either laying down or sitting upright with your eyes closed.

7. After use, ensure the power is off and the controller is disabled, wipe down both masks with a damp cloth and store flat.

8. For Tips & Tricks, visit puredailycare.com

BENEFITS

Designed to enhance the appearance of your skin by addressing your skin concerns with LED light, this portable device uses light therapy to rejuvenate both your face and neck while enhancing the results of your topical skincare. With seven different cool light functions to choose from, you can personalize your weekly treatments to what fits you best.

The Luma Mask Professional was inspired straight from Dermatology clinics and transformed into a world-renowned device that can be used whenever you want, from the comfort of your home. Lightweight and easily contoured to your face, the Luma Mask Professional can tackle all of your aging and skincare concerns.

Concern: Aging

Light: Red (630nm)

Benefit: Combats and fights wrinkles and fine lines. Red light also helps stimulate skin's collagen production and promotes circulation, for smoother and firmer skin.

Concern: Acne

Light: Blue (415nm)

Benefit: Fights acne-causing bacteria while reducing blemishes and aiding in oil control.

Concern: Pigmentation

Light: Green (525nm)

Benefit: Helps to even the skin tone and address pigmentation. Helps fight age spots dark patches and more.

Concern: Redness

Light: Yellow (570-590nm)

Benefit: Provides relief to skin's redness caused from sensitive skin, breakouts, dryness, flakiness, and overall inflammation.

Concern: Dull Lackluster Skin

Light: Purple (390-410nm)

Benefit: Promotes increased cell rejuvenation and revitalizes skin with a brighter and more radiant glow.

Concern: Irritation

Light: Cyan (480-490nm)

Benefit: Soothes irritated skin and helps reduce the look of capillaries and splotchy skin.

Concern: Overall Skin Benefits

Light: White (all)

Benefit: Amplifies skincare treatments by allowing products to penetrate deeper into the skin.

